


LEMARC

U S A

Model LK21 electric sonic toothbrush for children

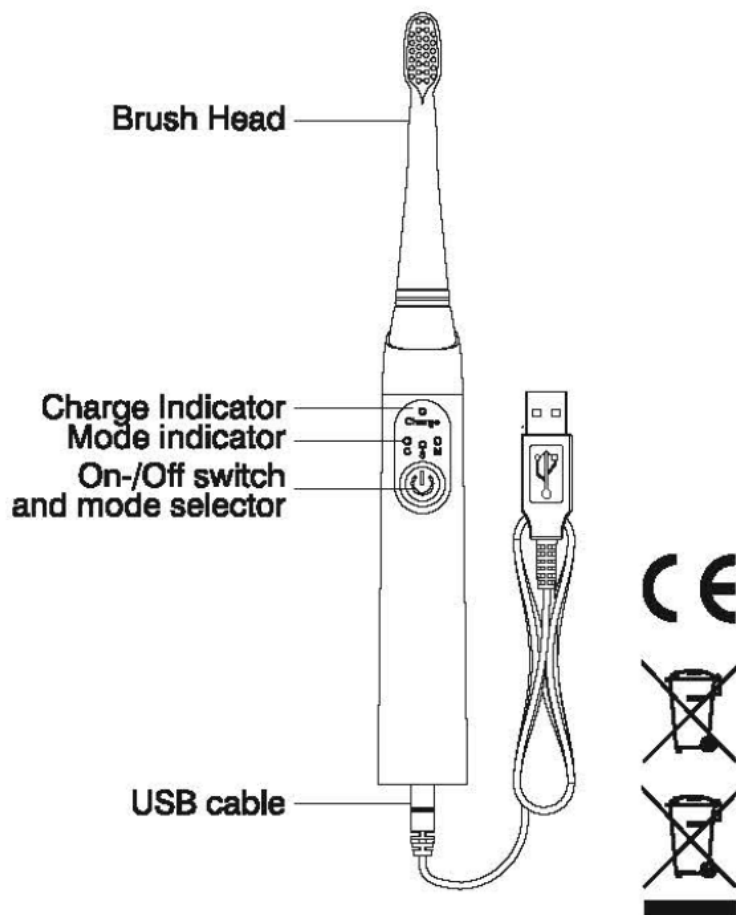


Charging


- All our electric toothbrush is delivered pre-charged, so you can use them right away. In case of longer storing, the battery might be exhausted. In this case charge first.
- Please charge your electric toothbrush by connecting it with the supplied cable to any USB port, for example the charger of your mobile phone or tablet, or on your computer.
- While charging the LED of the charging indicator blinks. Once it stops blinking and lights permanently, charging is complete. You can use your brush now for up to 40 days, without charging again (based on 2 cycles of 2 minutes per day).

- It is normal, that the housing of the handset gets slightly warm during charging, which is usually in a temperature range comfortable to touch. However, in case it turns hot and can't be touched comfortably, disconnect the power immediately and contact us via support@lemarcusa.com

Getting Started

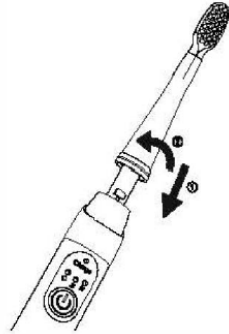


Please install the brush head as shown and turn it tight. Taking off the brush head follow these steps vice versa.

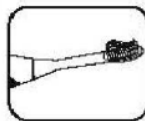

Moisten the bristles and apply the required amount of toothpaste.




Place brush head in the mouth before switching the appliance on. Place the bristles on your teeth in a 45° angle, facing the gum. Turn the brush on.

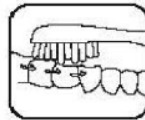

Now move the brush head with a wavy pattern along your teeth, applying a little pressure. While cleaning, the middle part of the bristles should always be in contact with teeth. Do not scrub, let the brush do the work for you.

Smart Timer

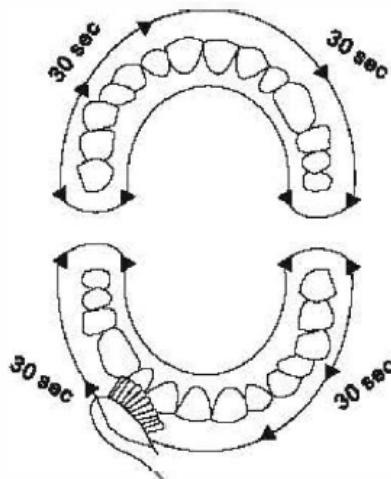

The toothbrush features a smart 2-minutes timer function. After each 30 seconds, a small interruption of the vibrations will remind you to move to the next section of your teeth. After 2 minutes the toothbrush will turn off automatically. 2 minutes is medically correct brushing, without watching the clock.

Functions

- This toothbrush features 3 modes **C** for children 7-11 years old , **S** (with reduced vibration) for kids 3-6, and **M** (with changing frequencies) for children 7+ to gums massaging.
- Pressing the On-/Off switch within 5 seconds after turning on, will change the mode. Pressing it later than 5 seconds after turning on, will turn the brush off.
- After switching on again, the appliance starts in the last selected mode.

Warnings

- This toothbrush is for personal use only. It is not suitable for commercial or professional use on multiple people, for example in dental Clinics.
- Use of the same brush head by several persons is not hygienic and may lead to transfer of infections and diseases. The brush head should be changed latest after each 3 months period.
- This toothbrush is ONLY suitable for cleaning teeth, gum and tongue. It must not be used for any other purpose and must not be inserted inside the body or any opening thereof. It is not a toy.
- Should using this toothbrush cause you any discomfort, pain, heavy bleeding or bleeding over a longer period, stop using the brush and contact your doctor.
- Within 2 months after a dental or gum surgery, please contact your doctor before use.
- This toothbrush has a reduced frequency and is suitable for children 3+ y.o. Children between 6 and 10 years should be instructed before or monitored while using. Children below 6 years must be monitored while using.
- The toothbrush must not be used with the charging cable connected to it.
- The toothbrush may not be used if brush head, handset or charging cable are damaged. The handset may only be opened and maintained by professional personal of authorised service suppliers. In case the battery needs to be changed or any other maintenance or repair works are required, please contact with our appointed service points.
- The operation temperature of this device is from 0° to 40°C. Please be reminded, that the battery exhausts significantly faster at lower temperatures.
- Clean the brush thoroughly after each use. Flush the brush head in water and clean the handset with a humid cloth.
- For a long time of your toothbrush, please store it in a dry place and do not expose it to direct sunlight for longer time. The brush must not come in contact with corrosive, abrasive, etching or solvating substances and materials, such as detergents, acids, solvents, etc.
- This product, as well the battery contained therein, must not be disposed off with normal household rubbish (2012/19/EU and 2006/66/EC). Follow your country's rules for the separate collection and recycling of electrical and electronic products and rechargeable batteries. Correct disposal helps to protect the environment.

Features



2 – Minutes Timer



Waterproof IPX7



3 modes



Charge once for up to 40 days



Charging via USB



Quality bristles by Dupont

For all
Service enquiries
Please contact us: support@lemarcusa.com

This a quality product by

LEMARC
U S A

WARRANTY CARD

Dear customer, thank you for purchasing LEMARC USA products!

After buying a product, please check whether the accessories are complete. Before turning on the device, please carefully read the operating instructions. The warranty period for LEMARC USA equipment is 12 months from the date of purchase. Corrections, deletions or other changes in the warranty card are not allowed and lead to loss of warranty service.

With this warranty card, the manufacturer confirms that it has assumed the obligation to satisfy the consumer's requirements established by the current consumer protection legislation in the event that defects in the product are discovered. At the same time, the manufacturer reserves the right to refuse warranty service of the product in case of non-observance of the conditions set forth above.



support@lemarcusa.com

www.lemarcusa.com